



箱型乾燥機

OVEN DRYER

+

+



原理介紹 +

將所要乾燥的產品,鋪放於一個盤子上,且將盤子放入可活動的台車,分成多層擺放,利用熱空氣將產品內水份帶走,而達到乾燥的目的。


適用產業 +

西藥、中藥、生物科技、健康食品、化工產業之各種濕顆粒、粉末之乾燥。

功能 +

將含有水份的粉末,顆粒,細顆粒,長條顆粒,小圓粒...等,乾燥至所需之含水量。

流程圖 | FLOWCHART +

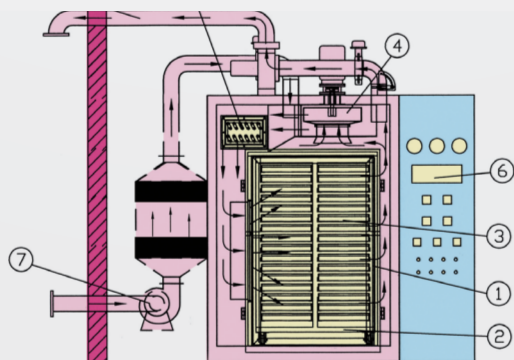

- | | |
|--------------------|--------------------|
| 1. Drying Chamber | 5. Heat Exchanger |
| 2. Drying Trolley | 6. Control Panel |
| 3. Drying Tray | 7. Fresh Air Inlet |
| 4. Circulation Fan | 8. Exhaust Damper |

THEORY +

Put materials on the trays with screens or perforations, and put the trays by layers on an active trolley. The hot air inside the dryer would remove moisture and dry the products.

APPLICATION +

Drying for pharmaceutical, Chinese herbal medicine, bio-tech product, Health food, any kind of chemical wet granule and powder.

FUNCTION +

Drying for powder, granule, fine granule, long granule, pellet which contains moisture until get required content of moisture.

規格表 | SPEC

型號 Model	產能/批次 Capacity/Batch (Kgs)(B/D:0.5)	尺寸 Size(WxLxH) (MM)
YS-DO-1	60-100	800X2200X1950
YS-DO-2	75-150	800X3000X1950
YS-DO-3	100-200	800X3200X2200
YS-DO-4	150-300	1200X3500X2500
YS-DO-5	200-400	1200X3500X2800